


Headwall’s Hyperspec® III software brings together a powerful set of hyperspectral acquisition and data manipulation tools in a single, easy-to-use environment for use across a wide range of applications.


Hyperspec III software is extremely versatile, using the industry-standard CameraLink and USB interfaces while supporting Windows 7, Windows 8, and Linux operating systems. It implements seamlessly into everything from advanced machine vision systems to airborne applications and interfaces with all of Headwall’s hyperspectral sensors. Across all operating situations, Hyperspec software provides high frame rates in excess of 400 fps and is multi-threaded for simultaneous processing and scalability. In a machine-vision application, for example, Hyperspec software integrates with downstream robotics via FTP and socket messaging. For airborne applications, orthorectification is crucial when capturing hyperspectral data. Hyperspec software interfaces with GPS devices to make this possible and to allow for precise calibration and start/stop data capture (which can also be triggered using a time-based method). Also, multiple sensors can be simultaneously controlled by the software. Regardless of the application, Hyperspec III software manages all camera parameters and motion-control features such as pan-and-tilt as well as Headwall’s Starter Kit with moving web.

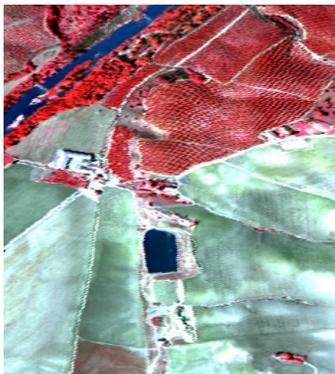

Airborne hyperspectral image

Because Hyperspec III software and all connected sensors can be controlled remotely, full automation is possible. Airborne remote sensing applications can be fully automated with respect to control and operation of the sensor. Working in conjunction with Hyperspec III is Headwall’s Hyperspec Data Processing Unit (HDPU), which is available in standard and lightweight airborne versions. Multiple display options exist within Hyperspec III software, including histogram and waterfall modes.

Application-Specific Solutions For Critical Environments

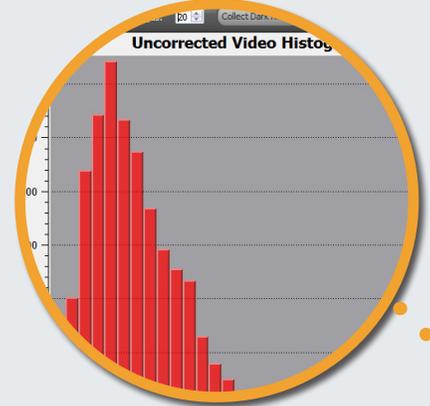

One software platform controls multiple sensors simultaneously!




Hyperspec[®] III Software


Headwall's powerful Hyperspec III software is designed for laboratory, ground-based, and airborne applications. The software is pre-loaded onto Headwall's Hyperspec Data Processing Units (standard and compact versions). Depending on the application and use-case, various modules within the software package are available. For example, airborne applications will benefit from GPS and Polygon perimeter triggers, plus ortho-rectification capabilities during post-processing.




Waterfall Display Settings

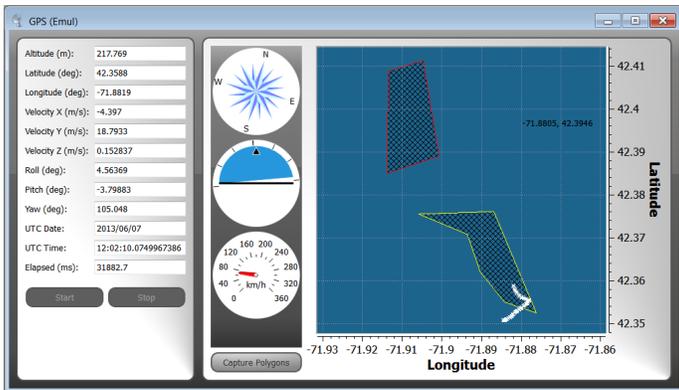


GPS

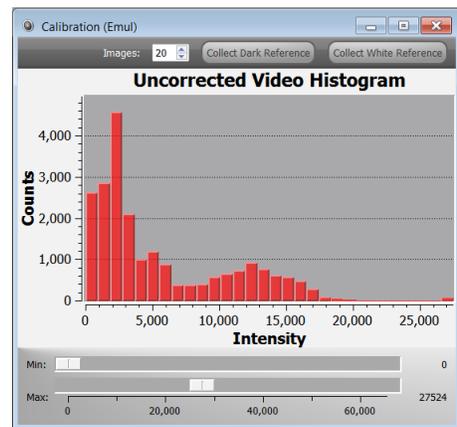

Histogram

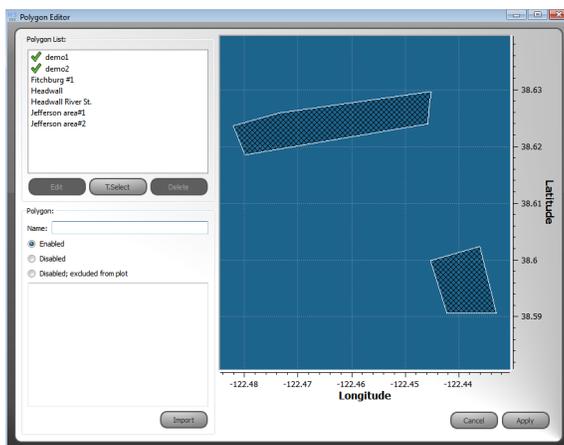

Polygon Editor

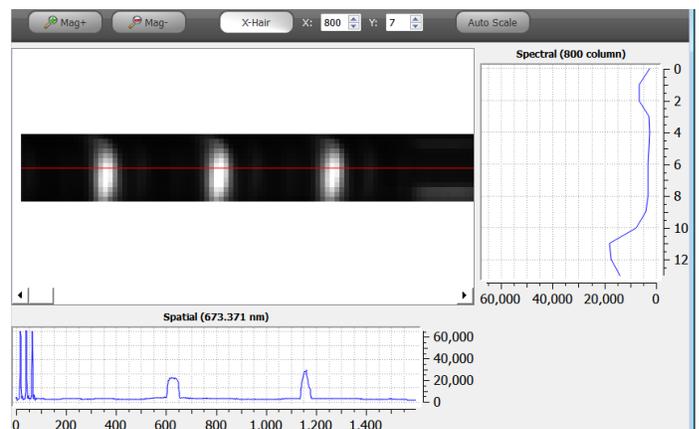

Live View

For more information contact:

analytikLtd (UK and Ireland Distributor)
 2 Cygnus Business Park, Middle Watch, Swavesey, Cambridge, CB24 4AA
 T: +44 (0)870 991 4044 E: info@analytik.co.uk www.analytik.co.uk



Headwall Photonics • 601 River Street • Fitchburg, MA 01420 • 978.353.4100 tel • www.headwallphotonics.com
 © Copyright 2013 Headwall Photonics, Inc. - The Hyperspec name is a registered trademark of Headwall Photonics, Inc.